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
```
cd /d %~dp0
```

```
call "C:\Program Files (x86)\LabSolutions\UVVis\Shimadzu\SpectroCommon\SystemBackup.exe"  
"C:\LabSolutions\LabSolutionsUVVis\UVVisSystemBackup.xml"  
pause
```

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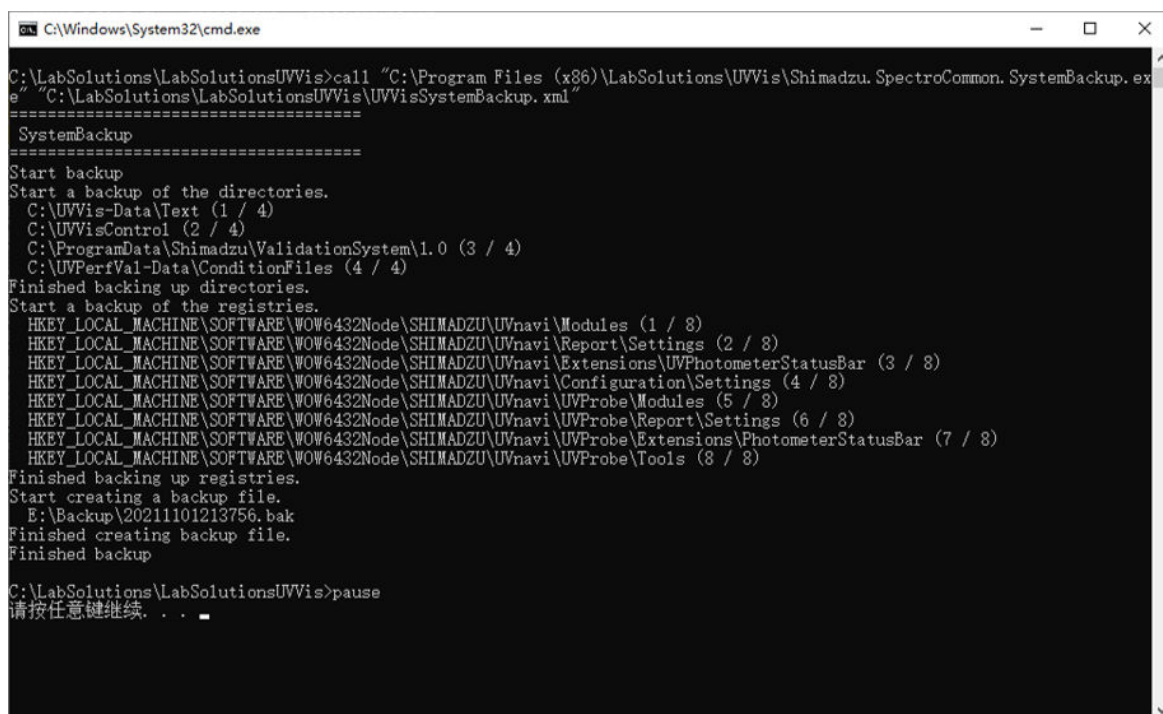```
C:\Windows\System32\cmd.exe

C:\LabSolutions\LabSolutionsUVVis>call "C:\Program Files (x86)\LabSolutions\UVVis\Shimadzu.SpectroCommon.SystemBackup.exe" "C:\LabSolutions\LabSolutionsUVVis\UVVisSystemBackup.xml"

=====
SystemBackup
=====
Start backup
Start a backup of the directories.
C:\UVVis-Data\Text (1 / 4)
C:\UVVisControl (2 / 4)
C:\ProgramData\Shimadzu\ValidationSystem\1.0 (3 / 4)
C:\UVPerfVal-Data\ConditionFiles (4 / 4)
Finished backing up directories.
Start a backup of the registries.
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Modules (1 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Report\Settings (2 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Extensions\UVPhotometerStatusBar (3 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Configuration\Settings (4 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Modules (5 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Report\Settings (6 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Extensions\PhotometerStatusBar (7 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Tools (8 / 8)
Finished backing up registries.
Start creating a backup file.
E:\Backup\20211101213756.bak
Finished creating backup file.
Finished backup

C:\LabSolutions\LabSolutionsUVVis>pause
请按任意键继续. . .
```

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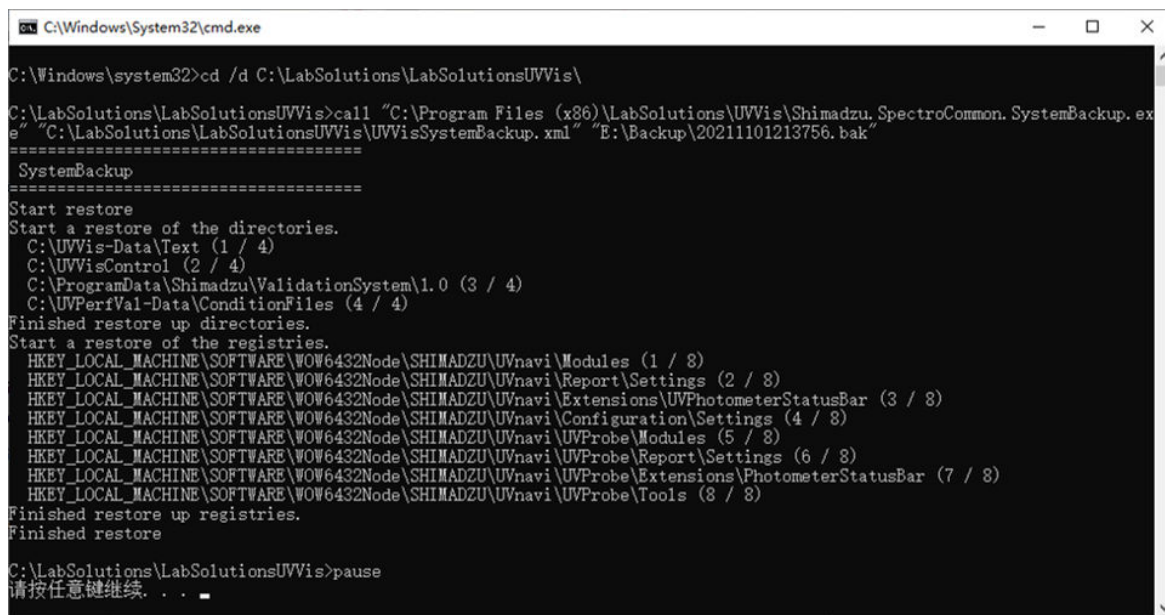
```
C:\Windows\system32>cd /d C:\LabSolutions\LabSolutionsUVVis\

C:\LabSolutions\LabSolutionsUVVis>call "C:\Program Files (x86)\LabSolutions\UVVis\Shimadzu.SpectroCommon.SystemBackup.exe" "C:\LabSolutions\LabSolutionsUVVis\UVVisSystemBackup.xml" "E:\Backup\20211101213756.bak"

=====
SystemBackup
=====
Start restore
Start a restore of the directories.
C:\UVVis-Data\Text (1 / 4)
C:\UVVisControl (2 / 4)
C:\ProgramData\Shimadzu\ValidationSystem\1.0 (3 / 4)
C:\UVPerfVal-Data\ConditionFiles (4 / 4)
Finished restore up directories.
Start a restore of the registries.
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Modules (1 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Report\Settings (2 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Extensions\UVPhotometerStatusBar (3 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\Configuration\Settings (4 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Modules (5 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Report\Settings (6 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Extensions\PhotometerStatusBar (7 / 8)
HKEY_LOCAL_MACHINE\SOFTWARE\WOW6432Node\SHIMADZU\UVnavi\UVProbe\Tools (8 / 8)
Finished restore up registries.
Finished restore
C:\LabSolutions\LabSolutionsUVVis>pause
请按任意键继续. . .
```

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